

A WORLD LEADER IN ORTHODONTIC APPLIANCES

OFFERS:

- Nite-Guide®
- Occlus-o-Guide®
- Interim-G®
- Ortho-T®
- Ortho-Tain Preformed Positioner®
- Habit-Corrector®
- Hawley-Form™
- Perfit Bumper®
- Snore-Cure®

APPLIANCES FOR ANY AGE TO CORRECT MOST MALOCCLUSIONS









AFTER







Nite-Guide®

Technique that guides incoming adult teeth in 93% of 5 to 7 year olds.

- Only TWO (2) appliances needed (one "C" Series and one "G" Series)
- Only ONE hour of passive nighttime wear needed, worn only while sleeping
- Guides incoming adult teeth into an ideal occlusion
- Scientifically designed sockets straighten erupting adult incisors





Guides in Adult Lower Centrals

- PDM fibers then form to hold aligned adult teeth straight permanently
- Improves the facial profile by mandibular advancement
- Corrects overjet, overbite, and many open-bites
- Orthodontic insurance covers this procedure
- "C" Series available in ELEVEN (11) sizes for optimum results

See page 15 for more details and measuring instructions.

EARLY TREATMENT FOR THE 5 TO 7 YEAR OLD





- **Severe Overjet Facial Disfiguration Sucked Thumb Severe Overbite**
- **ADULT DENTITION (12 YEARS, 8 MONTHS)**



Used Nite-Guide® ONLY while sleeping **Ideal Profile Ideal Occlusion** 5 years post-treatment

Potential Lower Crowding Extruded Upper Incisors Rotations





Straightened Adult Incisors Adequate Lower Space Good Facial Appearance









Excessive Open-Bite Sucked Thumb Speech Problems





Used Nite-Guide® ONLY while sleeping **Ideal Occlusion**

Occlus-o-Guide®

Mixed dentition treatment for 70% of malocclusions.

Corrects:

- Any severity of overjet
- Any severity of overbite
- Up to 4mm of crowding
- Up to 5mm of spacing
- Most TMJ problems
- Only one appliance needed for full treatment and retention
- Most cases finished in 4 to 10 months
- Takes only 1 to 2 hours of total chairtime
- Available in THIRTEEN (13) sizes for optimum fit



Occlus-o-Guide® appliance

See page 16 for more details and measuring instructions.

Occlus-o-Guide® ORTHO TECHNIQUE FOR THE 8-12 YEAR OLD








- **Severe Overjet Facial Disfiguration Class II Molars Deep Overbite**
- **ADULT DENTITION (14 YEARS, 4 MONTHS)**



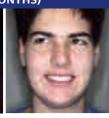

Used Occlus-o-Guide® 3 hours per day & nights **Ideal Occlusion Ideal Facial Appearance Retention result 5 years** after treatment

Crowded Upper Incisors Upper Incisors Rotated Facial Disfiguration





Used Occlus-o-Guide® 2+ hours per day & nights **Ideal Facial Profile Ideal Occlusion**





INITIAL DENTITION (9 YEARS, 11 MONTHS)





- Flared Upper Incisors **Severe Overjet Facial Disfiguration**
- **POST RETENTION (25 YEARS, 8 MONTHS)**





Occlus-o-Guide® worn 2+ hours per day & nights **Ideal Occlusion Ideal Facial Appearance**

Ortho-T®

Treats adults without the use of braces.

CORRECTS:

- Crowding
- Spacing
- Overjets
- Overhites
- TMJ

STRAIGHTENS ADULT CASES IN 1 TO 12 MONTHS

- Corrects rotations and crowding in 1 to 8 weeks
- Corrects up to 4mm overbite and overjet in 6 months

AFFORDABLE TREATMENT

- Only one appliance needed for treatment and retention
- No lab fees or waiting period prefabricated

DESIGNED TO ALSO CORRECT RELAPSED ORTHO CASES

WORN ONLY 1 TO 2 HOURS DAILY

- Wear it while driving to work, watching TV or surfing the web
- 22 hours a day with nothing on the teeth

- Only ONE removable appliance needed per patient
- No expensive lab fees
- Affordably priced treatment due to low overhead and minimal chairtime
- Correction in 1 to 12 months with as little as TWO (2) hours daily active wear
- Even severe crowding or spacing can be treated without stripping enamel when worn with a headgear or bumper
- THIRTEEN (13) preformed sizes available for optimal results.

See measuring instructions on page 17.

Ortho-T[®] ADULT BRACELESS ORTHO TECHNIQUE

















Ortho-T® CONTINUED

Ortho-T® can even correct severe cases when combined with a headgear or bumper with no stripping of precious adult enamel:

5mm of crowding • 5mm of overjet • 5mm of overbite









Impacted 2nd premolar
Deep overbite
Excess overjet
Ortho-T® worn with lower
bumper. Ortho-T® worn for
2+ hours per day and at night





Severe crowding with lack of space
Rotations
Ortho-T® worn only passively while sleeping
Cervical headgear worn at night simultaneously





Severe crowding
Rotations
Ortho-T® worn actively
2+ hours per day and at night
Bumper worn to distalize
lower molars







Crowding and rotations
End-to-end molar relation
Ortho-T® worn 2+ hours
per day actively
Headgear worn
at same time

Habit-Corrector®

TWO TYPES OF HABIT-CORRECTOR®

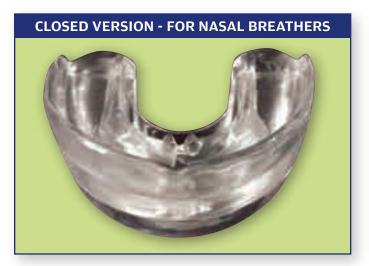

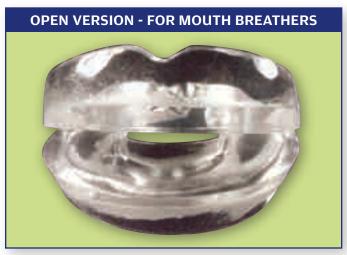

Aids in the correction of:

- Tongue Thrusting
- Abnormal Swallowing Patterns
- Thumb and Finger Sucking
- Open-Bites
- Overjets and Flared Incisors
- Mouth Breathing
- Speech Problems (Lisping)



Ortho-Tain[®] Positioners

ADVANTAGES:

- Immediate retention upon bracket removal
- Completely prefabricatedno lab work
- Lost appliances can be replaced immediately
- Retains overbite and overjet corrections
- Effectively corrects most TMJ problems
- Designed to provide orthodontic patients with an ideal Class Locclusion

- Also available in THREE (3) versions to accommodate: non-extraction cases
 (N series), (4) premolar extraction cases (X series), and (2) upper premolar extraction cases (U series)
- Available in various sizes.
- Rulers provided to determine correct size.
 For measuring instructions see page 17.

ORTHO-TAIN® appliances are made with Cooperation-Detector® (C-D®)

which turns the appliance slightly translucent when worn only at night, and becomes almost opaque white when worn 2 to 4 hours each day. If it is not worn for 48 hours, it returns back to transparent.

(Available for Occlus-o-Guide®, Ortho-T® and the Original Preformed Positioner®)



No Cooperation



Only Nighttime Wear



Daytime Exercise & Nighttime Wear

The Original Preformed Positioner®

TYPICAL ORTHO-TAIN® RESULTS













Interim-G Retainer®

- Used between Phase I & II of orthodontic treatment
- Worn only passively while sleeping
- Allows for deciduous exfoliation and adult eruption without adjustments
- Guides incoming adult teeth into perfect intercuspation

- 10 minute insertion– just measure and insert
- Extra soft material with breathing holes
- Available in clear and five colors: blue, magenta, green, purple and orange





INITIAL



PROGRESS - 1ST PHASE Interim-G Retainer® used in conjunction with a thumb-sucking appliance



END OF 1ST PHASE Interim-G Retainer® used while sleeping to retain and guide erupting teeth until beginning of 2nd phase ortho



FINAL Second phase of orthodontics finished

Hawley-Form™

THE 5-MINUTE CHAIRSIDE RETAINER



Appliances with Liner Kit

ADVANTAGES:

- Can immediately be customized in the mouth at chairside. No impressions or additional appointments needed
- Snaps firmly into place preventing appliance from falling out
- Maintains rotation corrections
- Retains overjets better than other retainers
- Can expand or constrict arch widths with the imbedded wire feature
- Can restore rotation relapses
- Can retain torque

The Hawley-Form™ appliances are available in five versions:







OR in combination with a lower for mandibular advancement:



Snore-Cure®



Snore-Cure® Appliance



Liner Kit

THE ADVANTAGES OF THE NEW SNORE-CURE® ARE:

- Can immediately be customized in the mouth at chairside. No impressions or additional appointments needed
- Snaps firmly into place to prevent falling out
- Advances mandible to prevent vibration of uvula
- Stays in all night
- One size fits all with hinge in back

The Softee® Appliance FOR THE 10 TO 12 YEAR OLD

Eliminates overbite and overjet before fixed appliance use Worn only while sleeping for about 4 to 6 months

- Worn only while sleeping to obtain overbite and overjet corrections
- Reduces fixed appliance treatment time up to one year
- Worn before and during bracket phase
- Obtains mandibular advancement for overjet correction
- Converts cases to Class I touch-ups using minimal brackets
- Super soft material for patient comfort
- Use as an anti-bruxism appliance
- Comes in 10 sizes

See page 16 for measuring instructions.




Perfit Bumpers®

FOR A "PERFECT" FIT

The Perfit® Bumper is a perfect adjunct to all our Ortho-Tain® appliances which are used for deciduous, mixed and adult dentitions that have 4 to 6mm of crowding. It is ideal for the adult, since the labial shield is hidden by the upper or lower lip, making it barely detectable.

- Distalizes deciduous or permanent molars (upper or lower) 1mm per month
- Takes only 1 to 3 try-ins for a perfect fit — saves
 15 minutes in fitting
- Comes in 8 convenient sizes and made of .045 wire
- Protects anchorage
- Protects mesial migration
- Labializes the lower or upper incisors if desired
- Can develop the lower and/or upper arch
- Plastic labial shield and hooks are adjustable with blue flame

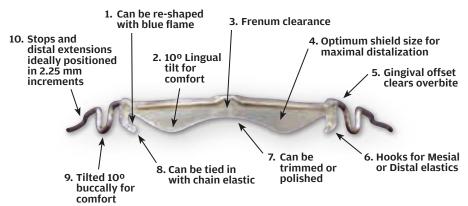


Sizes and Dental Arch Length from Molar Tube-to-Molar Tube (Mesial)

1L	57.0mm	5L	75.0mm
2L	61.5mm	6L	79.5mm
3L	66.0mm	7L	84.0mm
4L	70.5mm	8L	88.5mm

Most Common Range for Perfit® Bumpers

Adult Non-Extraction Sizes	2L - 7L
Mixed Dentition Sizes	3L - 8L
4 Bicuspid Extraction Sizes	1L - 6L
Deciduous Dentition Sizes	1L - 6L
Upper Arch Sizes	6L - 8L

Cleats for Rotation Corrections Available in TWO versions:

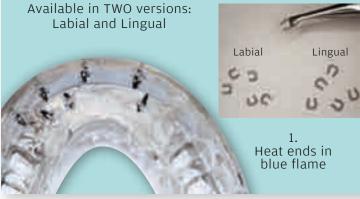





2. Insert into appliance

3.
Placed along long axis of tooth

Measuring Instructions

FOR THE 5 - 7 YEAR OLD

FOR USE WITH THESE APPLIANCES:

- Nite-Guide® "C" Series
- Nite-Guide® "G" Series (Occlus-o-Guide®)

The size of the appliance is stamped on upper left hand margin. Each appliance on the upper arch is 1.5mm apart and on the lower arch is 0.8mm apart.







To mesial of lower right canine



From mesial of upper left canine

HDDED



To mesial of upper right canine

Place pointer of pink ruler at mesial of lower or upper left deciduous canine; bend ruler along incisal edges (or along ideal arch form) and read size at mesial of right deciduous canine. Add two sizes to measurement for potentially crowded cases (example: $2 \rightarrow 2-1/2 \rightarrow 3$). This larger sized appliance is used in anticipation of expansion of the arches to provide space for the predicted adult crowding. The second appliance ("G" Series) is usually 3 sizes larger than the first ("C" Series) appliance used after the lower laterals begin erupting.

OR

"C" SERIES — FIRST APPLIANCE

First appliance ("C" Series) inserted when lower centrals begin erupting. Select "C" Series appliance two sizes larger than measurement (example: $2 \rightarrow 2-1/2 \rightarrow 3$).

"G" SERIES - SECOND APPLIANCE

Second appliance ("G" Series) inserted when lower laterals start to erupt. Select "G" Series appliance three sizes larger (example: $3C \Rightarrow 3-1/2 \Rightarrow 4 \Rightarrow 4-1/2G$).

Note: If the 6 year old molars are present when treatment is initiated, substitute a "G" Series appliance instead of a "C" Series appliance to prevent molars from super-erupting. The second appliance used when adult lower laterals start to erupt should be a "G" Series appliance that is 3 sizes larger than the first "G" Series appliance used. Both "G" Series appliances are worn only while sleeping.

"C" SERIES FOR DECIDUOUS NON-EXTRACTION CASES

	UPPER DECIDUOUS OR PERMANENT INCISAL AREA	LOWER DECIDUOUS OR PERMANENT INCISAL AREA	LOW CANINE TO CANINE WIDTH	%
SIZES	(MM)	(MM)	(MM)	POPULARITY
1 C	24.9-26.0	18.9-19.6	22.6	5%
1-1/2 C	26.1-27.2	19.7-20.4	23.2	10%
2 C	27.3-28.3	20.5-21.1	23.7	15%
2-1/2 C	28.4-29.5	21.2-21.9	24.3	15%
3 C	29.6-30.7	22.0-22.7	24.9	10%
3-1/2 C	30.8-31.9	22.8-23.5	25.4	8%
4 C	32.0-33.1	23.6-24.3	26.0	7%
4-1/2 C	33.2-34.2	24.4-25.0	26.6	6%
5 C	34.3-35.4	25.1-25.8	27.2	5%
5-1/2 C	35.5-36.6	25.9-26.6	27.7	3%
6 C	36.7-37.8	26.7-27.4	28.3	1%


(See "G" Series table on page 16)

Measuring Instructions

FOR THE 8 - 12 YEAR OLD

FOR USE WITH THESE APPLIANCES:

- Occlus-o-Guide® (G Series)
- Interim G® (IG Series)
- Softee® (CS Series)

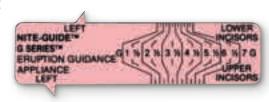



Place pointer at distal of upper left lateral



To distal of upper right lateral



From distal of lower left lateral



To distal of lower right lateral

Place the pointer of the pink ruler at the distal of the upper or lower left lateral and bend the ruler along the incisal edges to the distal of the right lateral. Be sure to use the correct side of the ruler marked either "lower incisors" or "upper incisors."

OR

Each size on the pink ruler is 1.5mm apart on the upper and 0.8mm on the lower. If there is no crowding or spacing the size measured on the ruler is used. If a patient has 1.5mm of upper anterior crowding and the size measured is a 4G, a 4-1/2G would be used to provide space for correction. If there is 3mm of upper anterior spacing and a 4G is measured on the ruler, a 3G would have to be used to close the spaces. In such a case, the appliance would be stretched across the spaced anterior teeth.

An alternate method is to measure the mesio-distal widths of the four upper or lower incisors in millimeters and select the Occlus-o-Guide® appliance that corresponds to this total measurement from the chart.

"G", "IG" or "CS" SERIES FOR MIXED-DENTITION CASES


	UPPER PERMANENT	LOWER PERMANENT	UPPER 6 PERMANENT	
SIZES	INCISORS (MM)	INCISORS (MM)	ANTERIORS (MM)	% POPULARITY
1 G	24.9-26.0	18.9-19.6	39.5-40.9	2%
1-1/2 G	26.1-27.2	19.7-20.4	41.0-42.4	3%
2 G	27.3-28.3	20.5-21.1	42.5-43.9	4%
2-1/2 G	28.4-29.5	21.2-21.9	44.0-45.4	6%
3 G	29.6-30.7	22.0-22.7	45.5-46.9	7%
3-1/2 G	30.8-31.9	22.8-23.5	47.0-48.4	10%
4 G	32.0-33.1	23.6-24.3	48.5-49.9	14%
4-1/2 G	33.2-34.2	24.4-25.0	50.0-51.4	16%
5 G	34.3-35.4	25.1-25.8	51.5-52.9	14%
5-1/2 G	35.5-36.6	25.9-26.6	53.0-54.4	12%
6 G	36.7-37.8	26.7-27.4	54.5-55.9	10%
6-1/2 G	37.9-39.0	27.5-28.2	58.0-57.4	3%
7 G	39.1 & over	28.3 & over	57.5 & over	2%

Measuring Instructions

FOR OVER 12 YEARS (ADULT)

FOR USE WITH THESE APPLIANCES:

- Positioner® (N Series Non Extraction)
- Positioner® (X Series 4 Premolar Extractions)
- Positioner® (U Series 2 Upper Premolar Extractions)
- Ortho-T®

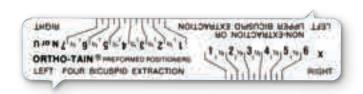





From distal of upper left canine

To distal of upper right canine

Place the pointer of the white ruler at the distal of the upper left canine then bend it along the incisal edges to the distal of the upper right canine. For retention after orthodontic bracket removal select an appliance that is one size smaller than the size indicated on the ruler for an optimal result.

Note: One side of the ruler measures non-extraction and two upper premolar extractions. The other side measures four premolar extractions.

"N" SERIES NON-EXTRACTION CASES

SIZES	UPPER 6 PERMANENT ANTERIORS (MM)	% POPULARITY
1 N	39.5-40.9	2%
1-1/2 N	41.0-42.4	3%
2 N	42.5-43.9	4%
2-1/2 N	44.0-45.4	6%
3 N	45.5-46.9	7%
3-1/2 N	47.0-48.4	10%
4 N	48.5-49.9	14%
4-1/2 N	50.0-51.4	16%
5 N	51.5-52.9	14%
5-1/2 N	53.0-54.4	12%
6 N	54.5-55.9	10%
6-1/2 N	56.0-57.4	3%
7 N	57.5 & over	2%

"X" SERIES FOUR (4) BICUSPID EXTRACTION CASES

SIZES	UPPER 6 PERMANENT ANTERIORS (MM)	% POPULARITY
1 X	45.5-46.9	3%
1-1/2 X	47.0-48.4	6%
2 X	48.5-49.9	7%
2-1/2 X	50.0-51.4	9%
3 X	51.5-52.9	11%
3-1/2 X	53.0-54.4	16%
4 X	54.5-55.9	17%
4-1/2 X	56.0-57.4	12%
5 X	57.5-58.9	9%
5-1/2 X	59.0-60.4	6%
6 X	60.5 & over	4%

"U" SERIES TWO (2) UPPER BICUSPID EXTRACTION CASES

SIZES	PERMANENT ANTERIORS (MM)	% POPULARITY
2 U	42.5-43.9	4%
2-1/2 U	44.0-45.4	6%
3 U	45.5-46.9	8%
3-1/2 U	47.0-48.4	11%
4 U	48.5-49.9	16%
4-1/2 U	50.0-51.4	18%
5 U	51.5-52.9	16%
5-1/2 U	53.0-54.4	12%
6 U	54.5 & over	9%

Now on DVD...

OCCLUS-o-GUIDE® &

NITE-GUIDE® TECHNIQUES

#1 OCCLUS-o-GUIDE® 3 hours on Overbite and Overjet correction

#2 OCCLUS-o-GUIDE® 3 hours on Rotations, Crowding and TMJ

#3 NITE-GUIDE® 3 hours on early intervention in the

4 to 7 year old to prevent adult malocclusions

Earn 6 hours credit toward certification by the Academy of General Dentistry for the review of each DVD set for a total of 18 hours of credit!

Courses are offered in different cities throughout the year. Check online at www.ortho-tain.com for the most up-to-date course schedule.



WAITING ROOM TAPE

OCCLUS-o-GUIDE® and/or NITE-GUIDE® Promotional Tape for Waiting Room Viewing Fifteen minute presentations with continuous viewing for 2-1/2 hours

PATIENT INSTRUCTIONAL TAPE

Ten minute chairside instructional video explaining use, care and wear for new Occlus-o-Guide® patients.

MONOGRAPHS

Over 60 complimentary monographs available online at www.ortho-tain.com INSIGHTS ON USE • SYNOPOSIS OF TREATMENT • FREQUENTLY ASKED QUESTIONS

WEBSITE PRESENTATIONS

Available at www.ortho-tain.com on the Nite-Guide®, Occlus-o-Guide® and the Ortho-T® techniques.

Over 2 MILLION used to date! Use the Original Ortho-Tain® Preformed Appliances DON'T ACCEPT SUBSTITUTES!











Occlus-o-Guide®

Mite-Guide®



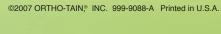

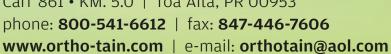

Carr 861 • KM. 5.0 | Toa Alta, PR 00953





AFTER















