
The Miniprevail Family Your artistry goes hand in hand with our design.



mini**Prevail** TWIN

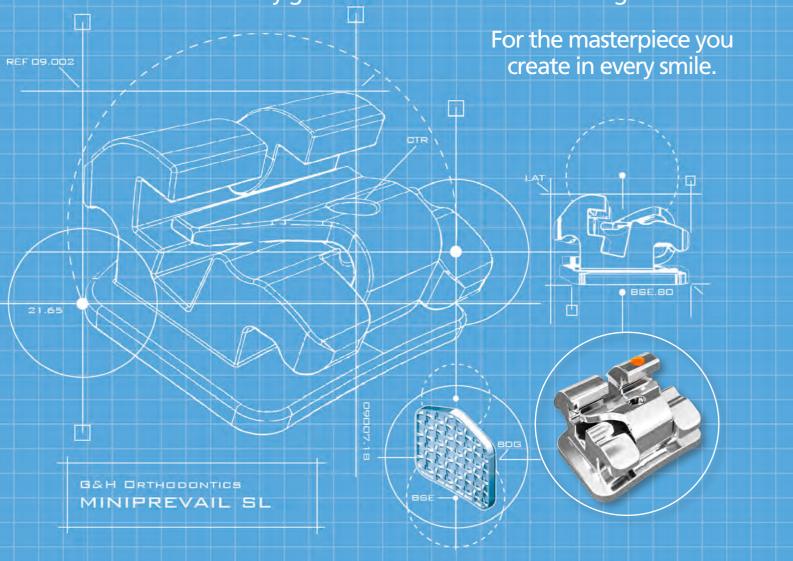

A low profile, miniature twin bracket that incorporates all the benefits of a larger bracket.

mini**Prevail TUBE**



Our distinct design takes the guesswork out of placement while the tube features offer easy wire placement and engagement.

mini**Prevail**'sL



An interactive self-ligating bracket. The answer to achieve excellent, predictable results in less time.


Shop our full line of G&H manufactured products @ **GHOrthodontics.com**





Our low profile bracket precision manufactured by G&H® Orthodontics


miniPrevail® TWIN brackets are constructed of a specialized micro-grain stainless steel MIM powder, ideal for creating small, robust appliances that are inherently stronger and more dimensionally accurate than other similar brackets.

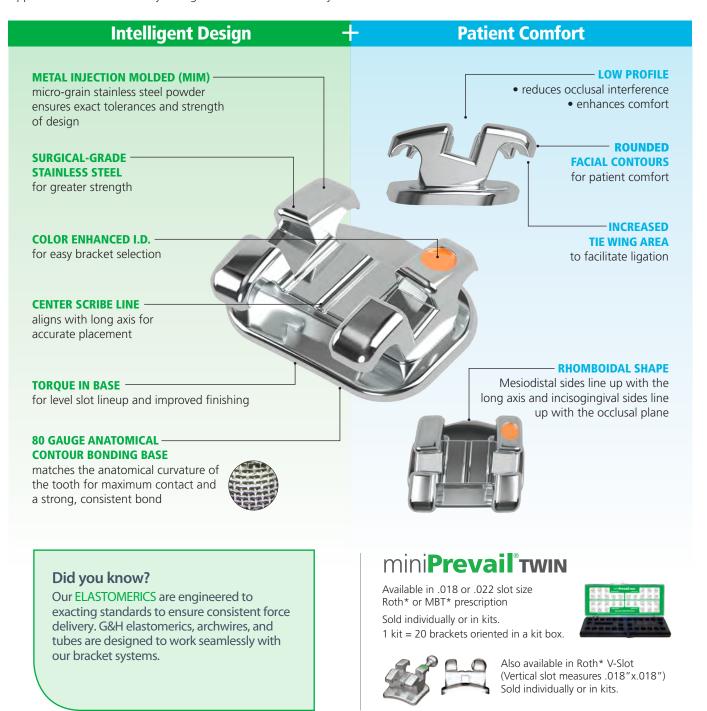



The answer to achieve excellent, predictable treatment results in less time.

An **INTERACTIVE** self ligating bracket which uses a NiTi spring clip that is strong enough to deliver continuous and consistent forces while assuring tooth movement control. When using round archwires (up to .020") allows **PASSIVE** engagement during initial leveling and alignment to reduce friction and improve sliding mechanics. When using rectangular archwires (from .016"x.022"), becomes **ACTIVE** to provide effective torque and rotation control toward the end of treatment.

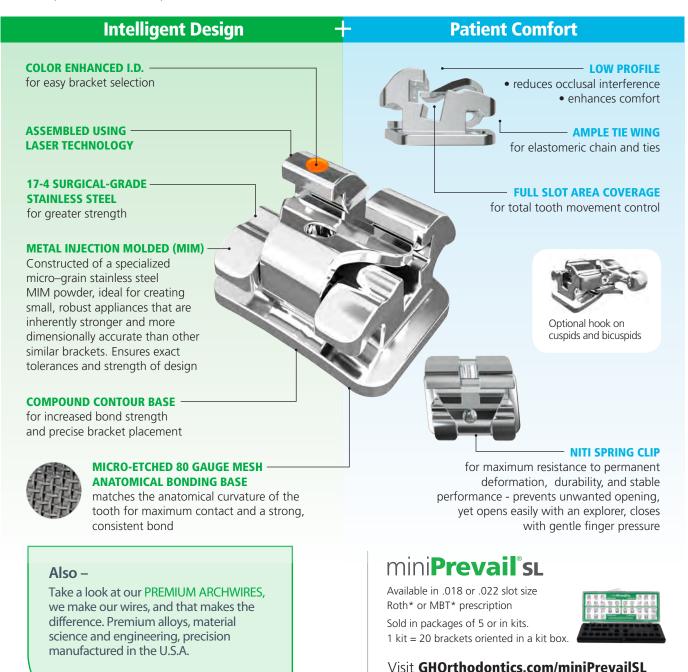


Shop our full line of G&H manufactured products @ **GHOrthodontics.com**

See our website or catalog for full details.



Seamlessly working together with our miniPrevail®TWIN and SL

Our new miniPrevail®TUBE has been engineered with the patient and orthodontist in mind. Our distinct "crosshair" design takes the guesswork out of placement while the tube features offer easy wire placement and engagement. Finally, the low profile and sleek hook is designed for the utmost patient comfort.

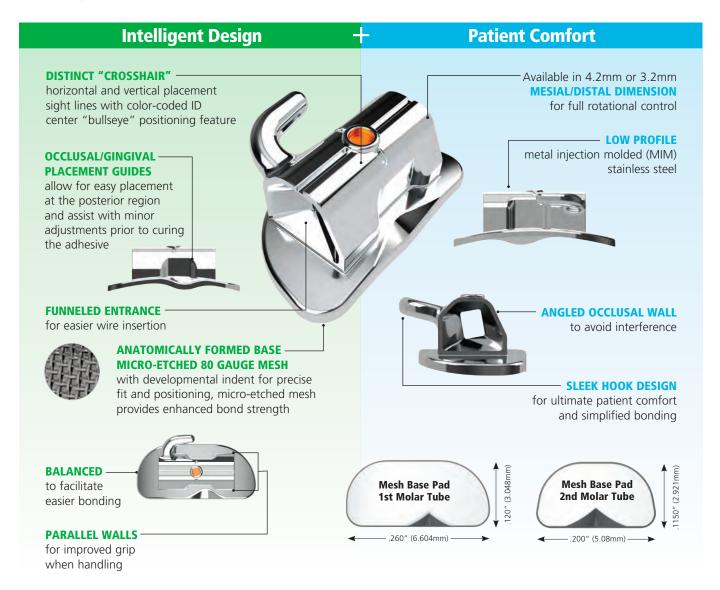

Need BANDS?

Our BANDS, provide a secure fit you can count on. Precision-shaped tooth anatomy for a tight, secure fit and excellent retention, time after time.

See our website or catalog for full details.

mini**Prevail** TUBE

Available in .018 or .022 Weldable or Bondable Roth*, MBT* or Standard prescription Sold in packages of 10.



