


Bonds Hooks and Buttons to Clear Aligners

NO SURFACE CONDITIONING REQUIRED.


- Will bond any button or hook to clear aligner materials - Essix, Invisalign, Invisalign Smart Track Bioplast, Tru-Tain, and more!
- A clear, light-cure paste for excellent aesthetics.
- Bonds to aligners without plastic preparation, primers, or conditioners.
- Low modulus of elasticity - Bond Aligner® flexes with thermoformed plastic trays and will not crack or separate.

ONLY FROM  **RELIANCE®**


Simple Application

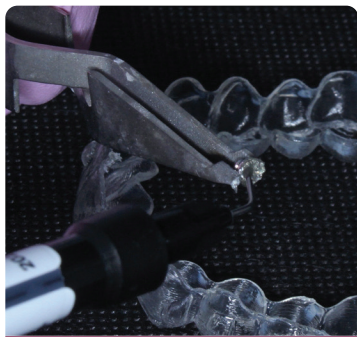

1. Apply Bond Aligner to bonding base.

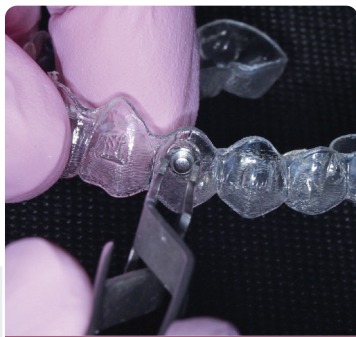

2. Deliver and position.

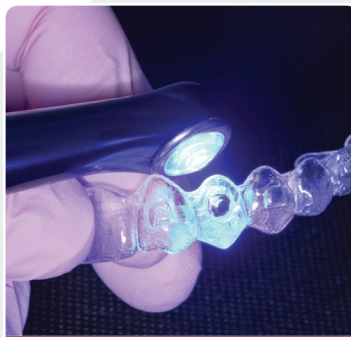

3. Light cure.



Bonded aligner.



Secure bonds, active application.



Reliance Orthodontic Products:
Over 40 Years of Serving the
Orthodontic Profession.



RELIANCE®
ORTHODONTIC
PRODUCTS

O (800) 323-4348
M (630) 773-4009

F (630) 250-7704

RelianceOrthodontics.com
RelianceU.tech

Proud Sponsor of
AAOF